
## Introduction to Project 5: MRI Measurements of Perfusion

Quantitative and Functional Imaging
BME 4420/7450
Fall 2022


# Review: magnetic susceptibility contrast agents

- Large magnetic moment
  - Gadolinium (7 unpaired electrons)
- Increase relaxation rate of neighboring water molecules
  - T<sub>1</sub> relaxation (M<sub>z</sub>)
  - $-T_2$  relaxation  $(M_{xy})$
- When confined to vessels, concentrated magnetic moments shift the field around the vessels

#### Magnetic field surrounding a vessel

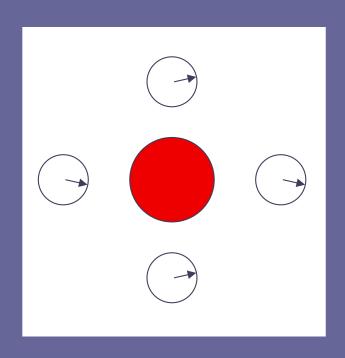


Normal blood

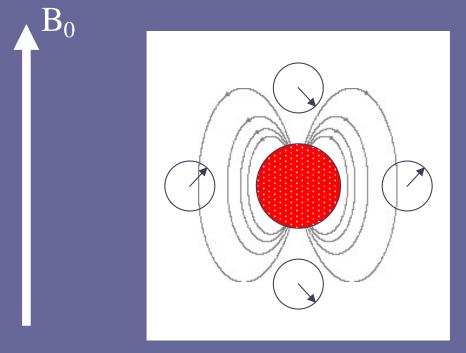

Blood with contrast agent

## Capillary networks

All cells are close to capillary network - fell magnetic effect of Ga contrast

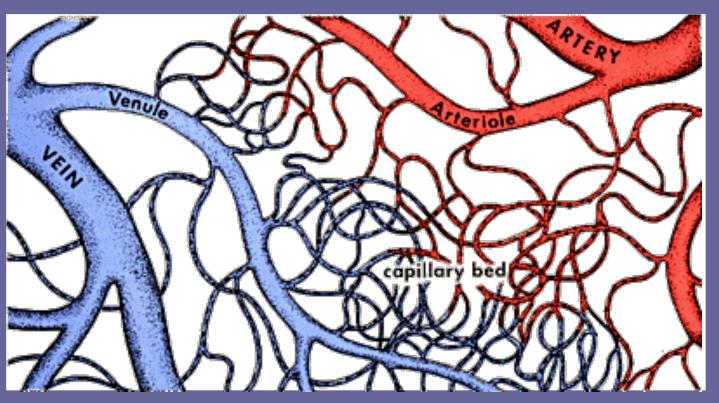

Buxton, 2002

## Dependence of relaxation rate on tissue contrast agent concentration

• Relaxation rate is the inverse of the time constant

$$R_{2}^{*} = \frac{1}{T_{2}^{*}}$$

$$R_{2}^{*} \text{ (with contrast)} = R_{2}^{*} \text{ (without contrast)} + \Delta R_{2}^{*}$$

• Tissue relaxation due to contrast agent

$$\Delta R_2^* = k \cdot C_T(t)$$

• Signal intensity prior to contrast agent injection

$$S_0 = A \cdot e^{-T_E \cdot R_2^* \text{(no contrast)}}$$

• Signal intensity with contrast agent


$$S = A \cdot e^{-T_E \cdot R_2^* (\text{contrast})}$$

$$= A \cdot e^{-T_E \cdot R_2^* (\text{no contrast})} \cdot e^{-T_E \cdot \Delta R_2^*}$$

$$= S_0 \cdot e^{-T_E \cdot \Delta R_2^*}$$

## Image analysis

- Find signal time course in a region of interest (ROI)
- Identify the baseline (pre-contrast R<sub>2</sub>\*)
- Convert signal changes to R<sub>2</sub>\* changes

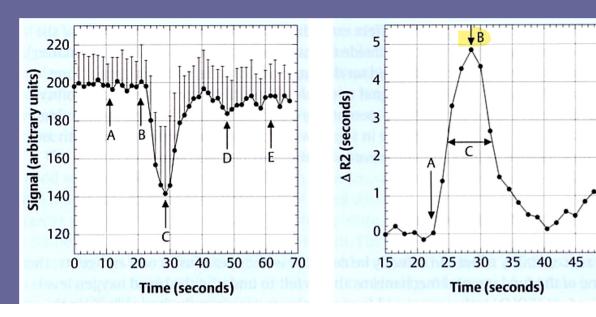

small peak = blood clot --large peak = normal

$$\Delta R_2^* = -\frac{\ln(S/S_0)}{TE}$$

Sorensen, 2000

### Image analysis

- Find signal time course in a region of interest (ROI)
- Convert signal changes to R<sub>2</sub>\* changes
- Identify the baseline (pre-contrast R<sub>2</sub>\*)
- Integrate the changes in R<sub>2</sub>\*

$$CBV = \frac{\int_0^\infty C_T(t)dt}{\int_0^\infty C_A(t)dt} = \mathbf{C} \cdot \int \Delta R_2^*(t)dt$$

### Project 5 Goals

- Measure image intensity changes in an ischemic region of interest (ROI) in the brain of a stroke patient
- Estimate the mean cerebral blood volume (relative CBV= rCBV) in the ROI

$$\Delta R_2^* = -\frac{\ln(S/S_0)}{TE}$$

$$rCBV = \int \Delta R_2^*(t)dt$$

– Is this significantly different from healthy tissue?

#### Problem #1: data visualization

- How good are the data?
  - Contrast
  - Noise
  - Artifacts
- Check the image time course visually
- Show images in a movie


## Making movies



#### Problem #2: what is a normal rCBV?

- We need to compare rCBV in the stroke region to a normal rCBV
  - Where can we get this?

right side signal reduce faster & recover faster ---> left side has sth that slow flow of blood BLOOD clot / stroke

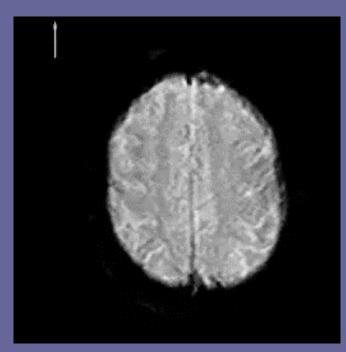

Sorensen, 2000

#### Problem #2: what is a normal rCBV?

- We need to compare rCBV in the stroke region to a normal rCBV
  - Where can we get this?
- Compare to a similar region in the opposite hemisphere